


Getting Started Guide

TD I-110




What's in the Box?

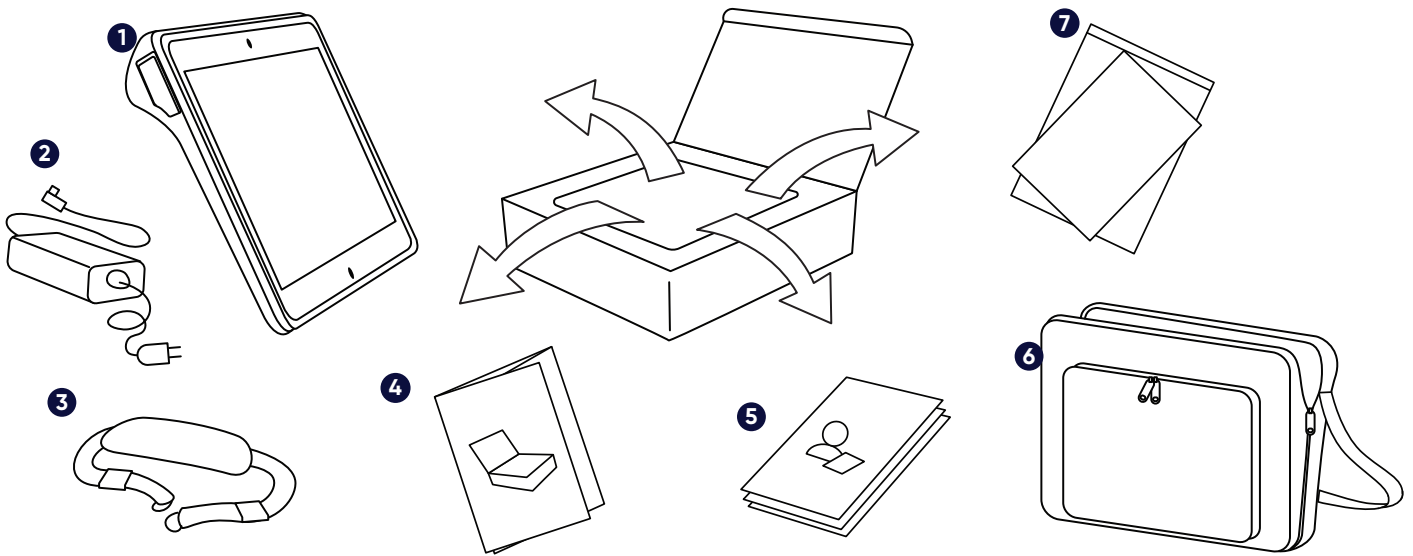

- 1 Device
- 2 Power Cable
- 3 Shoulder Strap
- 4 Getting Started Guide
- 5 TD Snap Training Cards
- 6 Carry Case
- 7 Safety & Compliance, License, and Configuration Documents

Getting to Know Your Device

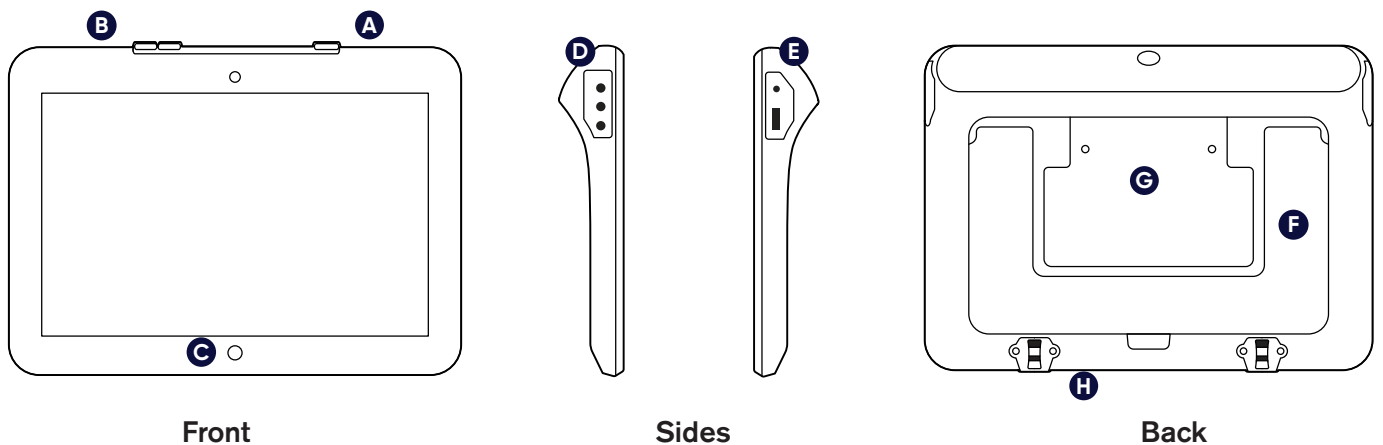



- A Power Button
- B Volume Buttons
- C Home Button
- D Audio & Switch Ports
- E Power Connector & USB Ports
- F Stand
- G Mounting Location
- H Strap Connectors

Follow the steps in this guide to set up and start using your device.

Step 1: Initial Setup

Start the Device


- 1 Connect the power cable to the device and plug in to an outlet
- 2 Press the Power button

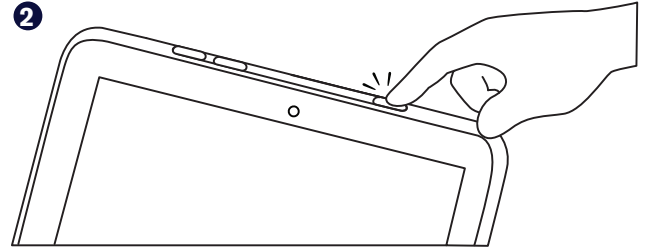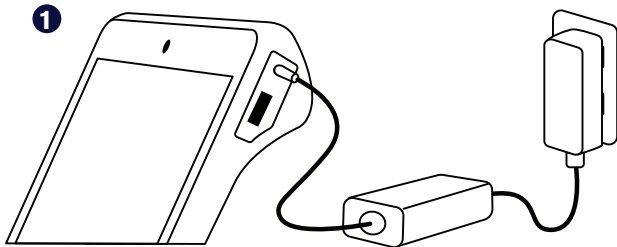

Windows Setup

If your device has not been pre-configured, you will be prompted to complete the Windows Setup process first.

While this process does not require a Microsoft account, we recommend that you connect to or create the user's account and not a caregiver's account.

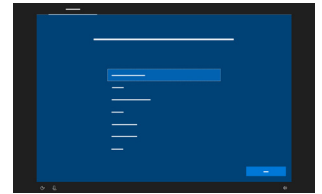

Step 2: Set Up TD Snap



Double tap the icon to launch the TD Snap app. Follow the prompts to create a new user or restore a saved user file.

Step 3: Mount & Position

The device may be positioned using a mounting system, laying it flat on a surface, or propped on the integrated kickstand. Start by situating the user comfortably, then find the device position that gives them clear screen visibility and easy access to their selection method of choice. Always position the device to suit the user, not the other way around. It is expected that the device may need to be repositioned throughout the day. Attach the shoulder strap for secure transport.

Step 4: TD Snap Training Cards

Continue your set up journey with the TD Snap Training Cards that came in the box with your device. The training cards walk you through the features of TD Snap, basic editing, backing up your data, and some activities to help you integrate TD Snap into your daily life.

Additional Resources

Scan the QR codes with your phone or use the links.



Tobii Dynavox

uk.tobiidynavox.com



TD Snap Help

qrco.de/bclHVF



**TD Facebook
Community**

qrco.de/TDFB



**Tobii Dynavox
Learning Hub (English
only)**

learn.tobiidynavox.com



**TD I-110 Device
Documentation**

qrco.de/I-110manual



myTobiiDynavox

mytobiidynavox.com



**Learn about access
methods**

qrco.de/bbA7US



UK Technical Support

0114 481 0011

support.uk@tobiidynavox.com